Panasonic®

Operating Instructions (Household use) Oral Irrigator Model No. EW1213

English.....	3
Español.....	19

Operating Instructions
(Household use) Oral Irrigator
Model No. EW1213

Contents

Safety precautions	5	Frequently asked questions	14
Important information	8	Troubleshooting	14
Parts identification	9	Battery life	15
Charging	9	Removing the built-in rechargeable battery	15
How to use	10	Specifications	16
Tips on usage	12	Changing the nozzle (Sold separately)	16
After use	12		
Cleaning	13		

Thank you for purchasing this Panasonic product.

Before operating this unit, please read these instructions completely and save them for future use.

IMPORTANT SAFEGUARDS

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following: **READ ALL INSTRUCTIONS BEFORE USING.**

DANGER To reduce the risk of electrocution:

1. Always unplug this product immediately after using.
2. Do not use while bathing.
3. Do not place or store product where it can fall or be pulled into a tub or sink.
4. Do not place in or drop into water or other liquid.
5. Do not reach for a product that has fallen into water. Unplug immediately.

WARNING To reduce the risk of burns, electrocution, fire or injury to persons:

1. Close supervision is necessary when this product is used by, on, or near children or invalids.
2. Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
3. Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return the product to a service center for examination and repair.

4

4. Keep the cord away from heated surfaces.
5. Never block the air openings of the product or place it on a soft surface, such as a bed or couch, where the air openings may be blocked. Keep the air openings free of lint, hair, and the like.
6. Never use while sleeping or drowsy.
7. Never drop or insert any object into any opening or hose.
8. Do not use outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.
9. Fill water tank with water only unless otherwise specified by manufacturer.

SAVE THESE INSTRUCTIONS

Safety precautions

To reduce the risk of injury, loss of life, electric shock, fire, and damage to property, always observe the following safety precautions.

Explanation of symbols

The following symbols are used to classify and describe the level of hazard, injury, and property damage caused when the denotation is disregarded and improper use is performed.

 DANGER	Denotes a potential hazard that will result in serious injury or death.
 WARNING	Denotes a potential hazard that could result in serious injury or death.
 CAUTION	Denotes a hazard that could result in minor injury or property damage.

The following symbols are used to classify and describe the type of instructions to be observed.

	This symbol is used to alert users to a specific operating procedure that must not be performed.
	This symbol is used to alert users to a specific operating procedure that must be followed in order to operate the unit safely.

WARNING

► Preventing accidents

-  **Do not store within the reach of children or infants. Do not let them use it.**
- Doing so may cause an accident or injury due to accidental ingestion of removable parts.

► Power supply

-  **Do not connect or disconnect the power plug to a household outlet with a wet hand.**
- Doing so may cause electric shock or injury.
- Do not immerse the charger in water or wash it with water.**
Do not place the charger over or near water filled sink or bathtub.
Do not use when the cord or the power plug is damaged or when the fitting into the household outlet is loose.
-  **Do not damage the cord or the power plug.**
 - **Do not deface, modify, forcefully bend or pull the cord or the power plug. Also, do not place anything heavy on, twist or pinch the cord.**- Doing so may cause electric shock or fire due to a short circuit.
-  **Do not use anything other than the supplied charger for any purpose. Also, do not use any other product with the supplied charger. (See page 9.)**
- Doing so may cause burn or fire due to a short circuit.

WARNING

Do not use in a way exceeding the rating of the household outlet or the wiring.

-  - Exceeding the rating by connecting too many power plugs to one household outlet may cause fire due to overheating.

Always ensure the appliance is operated on an electric power source matched to the rated voltage indicated on the charger.

-  **Fully insert the power plug.**
- Failure to do so may cause fire or electric shock.

► In case of an abnormality or malfunction

Immediately stop using and remove the power plug if there is an abnormality or malfunction.

- Failure to do so may cause fire, electric shock, or injury.

<Abnormality or malfunction cases>

-  • **The main unit or charger is deformed or abnormally hot.**
- **The main unit or charger smells of burning.**
- **There is abnormal sound during use or charging of the main unit or charger.**
- Immediately request inspection or repair at an authorized service center.

► This product

This product has a built-in rechargeable battery. Do not throw into flame or heat.

-  - Doing so may cause fluid leak, overheating, or explosion.

Do not modify or repair.

-  - Doing so may cause fire, electric shock, or injury.
Contact an authorized service center for repair (replacement of battery etc).

Never disassemble except when disposing of the product.

-  - Doing so may cause fire, electric shock, or injury.

► Cleaning

 **Always unplug the power plug from a household outlet when cleaning the charger.**

- Failure to do so may cause electric shock or injury.

Regularly clean the power plug to prevent dust from accumulating (once every half a year).

-  - Failure to do so may cause fire due to insulation failure caused by humidity.

Disconnect the power plug and wipe with a dry cloth.

► Note the following precautions

This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless

-  **they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.**

Children should be supervised to ensure that they do not play with the appliance.

- Failure to do so may cause an accident or injury.

-  **The supply cord cannot be replaced. If the cord is damaged, the charger should be scrapped.**

- Failure to do so may cause an accident or injury.

CAUTION

► Preventing damage to teeth and gums

Do not allow those that cannot operate the appliance to use this appliance.

 **People with weak oral senses should not use this appliance.**

-Doing so may cause injury or damage to teeth and gums.

People who may have periodontal disease, have teeth being treated, or concerned about symptoms within their mouth should consult a dentist before use.

-Failure to do so may cause injury or damage to teeth and gums.

► This product

Never use for a purpose other than oral cleaning.

-Doing so may cause an accident or injury when the stream of water is pointed toward the eyes, nose, ears, or throat.

Do not fill the tank with water warmer than 40 °C (104 °F).

-Doing so may cause burns.

Do not share nozzle with your family or other people.

-Doing so may cause infection or inflammation.

Do not place a metallic object such as a coin or a clip on top of the charger.

-Doing so may cause burns due to excessive heat.

► Note the following precautions

Do not allow metal objects or trash to adhere to the power plug.

-Doing so may cause electric shock or fire due to a short circuit.

Do not drop or subject to shock.

-Doing so may cause injury.

Do not wrap the cord around the charger when storing.

-Doing so may cause the wire in the cord to break with the load, and may cause fire due to a short circuit.

Disconnect the power plug from the household outlet when not using for extended period of time.

-Failure to do so may cause electric shock or fire due to electrical leakage resulting from insulation deterioration.

Disconnect the power plug by holding onto it instead of the cord.

-Failure to do so may cause electric shock or fire due to a short circuit.

After use, always care for and clean the appliance. (See page 13.)

-Failure to do so may cause damage to your health due to propagation of mold and bacteria.

Handling of the removed battery when disposing

English

DANGER

- The rechargeable battery is exclusively for use with this appliance. Do not use the battery with other products.
- Do not charge the battery after it has been removed from the product.
- Do not throw into fire or apply heat.
- Do not solder, disassemble, or modify the battery.
-  Do not let the positive and negative terminals of the battery get in contact with each other through metallic objects.
- Do not carry or store the battery together with metallic jewellery such as necklaces and hairpins.
- Never peel off the tube.
 - Doing so may cause fluid leak, overheating, or explosion.

The battery contains alkaline fluid. If it comes in contact with the eyes, do not rub the eyes, and thoroughly rinse with clean water, such as tap water.

- Failure to do so may cause loss of vision. Consult a physician immediately.

WARNING

After removing the rechargeable battery, do not keep it within reach of children and infants.

-  - The battery will harm the body if accidentally swallowed.
If this happens, consult a physician immediately.


The battery contains alkaline fluid. If it comes in contact with the skin or clothes, rinse off with clean water, such as tap water.

- Failure to do so may cause skin injury.

Important information

- Do not put salt, mouthwash or other chemicals into the water tank. Doing so may cause malfunction.
- Do not dry the main unit with a dryer or a fan heater. Doing so may cause malfunction or deformation of the parts.
- Do not wipe with thinners, benzine, alcohol etc. Doing so may cause malfunction, cracking or discoloration of the parts.
- Store the appliance in a place with low humidity after use. Leaving the appliance in a bathroom may cause malfunction.
- Charge the appliance correctly according to these Operating Instructions. (See page 9 “Charging”.)
- This appliance contains battery that is only replaceable by skilled persons. Please contact an authorized service center for the details of repair.

Parts identification

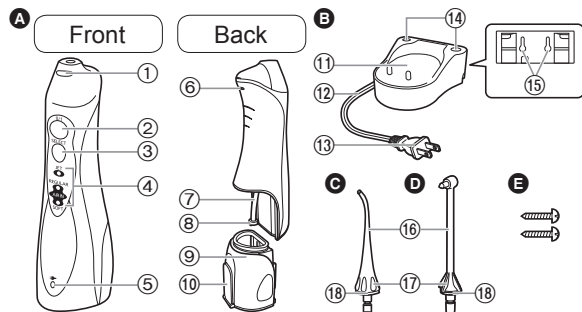A Main unit

- ① Nozzle release button
- ② Power switch [0/1]
- ③ Mode switch
- ④ Mode indicator
- ⑤ Charging indicator
- ⑥ Ventilation hole
- ⑦ Suction hose
- ⑧ Filter
- ⑨ Water tank
- ⑩ Water tank cap

B Charger (RE8-53)

- ⑪ Charging section
- ⑫ Cord
- ⑬ Power plug

C Nozzle (for general use)

- ⑭ Nozzle stands
 - ⑮ Screw holes for wall mounting
- The charger can be mounted on a wall using two wood screws.

D Orthodontic nozzle

- ⑯ Nozzle
 - ⑰ Nozzle handle
 - ⑱ Identification ring
- E 2 screws**


Note

There may be moisture inside the main unit (water tank and suction hose) left from the distilled water used for product testing.

Charging

► On the bathroom countertop or a shelf


1. Place the cord in the right, left or center cord guide, depending on where the charger will be placed.
2. Place the charger on a secure, level surface.
3. Plug in the power plug.

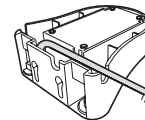1 Insert main unit so it stands upright in the charging section.

The charging indicator will light up.

- Be sure to place the main unit upright. If the main unit is tilted or not in direct contact with the charger, the main unit may not charge properly.
- Charging is completed in about 15 hours.
- The lamp does not go off even when charging ends.
- The main unit will provide an operation time of approximately 15 minutes of use when the battery is fully charged. (The operating duration will become shorter as the batteries get older.)

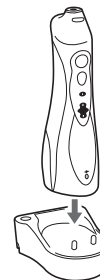2 Disconnect the power plug after charging is completed.

(for safety and to reduce energy usage)

- The battery performance will not be affected even if it is charged for more than 15 hours.

Notes

- When using the appliance for the first time or if it has not been used for more than half a year, the charging indicator may not light up for several minutes after charging starts or the operating time may become shorter. In such cases, charge the battery for 23 hours or more.
- Recommended ambient temperature for charging is 0 °C – 35 °C (32 °F – 95 °F). Battery performance may decrease outside of the recommended temperature.

How to use

1 Attach the nozzle firmly onto the main unit.

2 Open the water tank cap, hold the main unit so it is level, and fill the tank with water.

- Always use new tap water or lukewarm water.
- If the water feels cold, use lukewarm water with a temperature of 40 °C (104 °F) or less.

3 Close the water tank cap, making sure it locks in place.

4 Press the mode switch and select the desired water jet mode.

- The mode indicator will glow the previously selected mode. If the appliance is not used for more than 30 days, or after the appliance has been fully charged from its fully discharged state, the mode is initialized to AIR IN (SOFT). (When the battery is drained, the mode will be reset but this is not a fault.)
- If you want to use a different mode, press the mode switch until it changes to the desired setting. Pressing the mode switch will change the mode from AIR IN (SOFT) to AIR IN (REGULAR) to JET.

Mode	Description
	<p>JET mode</p> <ul style="list-style-type: none"> This mode is for removing food particles trapped between the teeth.
	<p>AIR IN (REGULAR) mode</p> <ul style="list-style-type: none"> This mode is for periodontal pocket cleaning and gum care.
	<p>AIR IN (SOFT) mode</p> <ul style="list-style-type: none"> This mode is for gently massaging the gums.

Notes


- When you use the appliance for the first time, press the power switch to turn the appliance on and test each of the water jet modes with the tank filled with water.
- The appliance may not operate in an ambient temperature approximately 5 °C (41 °F) or lower.
- Do not hold the tip of the nozzle firmly against the teeth or gums. (Except orthodontic nozzle)
- Since plaque adhering to the surface of teeth cannot be removed only with water flow, brushing is always necessary. It is recommended to clean inside your mouth with water flow after brushing your teeth.
- Do not operate the main unit when empty other than when pouring out any remaining water after use. Operating the main unit without water may result in a malfunction.

5 Place the nozzle in your mouth and close your mouth gently before pressing the power switch.

- To use, place your face over a sink, and with the main unit upright and the nozzle at the top, raise your elbow.
- Pressing the mode switch during use will change the mode.
- A fully filled water tank provides approximately 35 seconds of use. (When using the nozzle (for general use))
- If the nozzle is inclined too much, water will come out of the ventilation hole instead of the nozzle. (The ventilation hole in the rear of the main unit allows for easier water suction of the hose.)

Turn the nozzle handle to adjust the direction of the water flow.

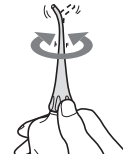6 After use, press the power switch to turn off the main unit.

- Stop operation with the nozzle still in the mouth to prevent water from splashing.

Tips on usage

Nozzle (for general use)

When cleaning with water flow

► Cleaning the space between the teeth

In the JET mode

- Direct the water flow onto the Δ zone between the teeth.
- Direct from the front and rear sides of the teeth.

► Periodontal pocket cleaning

In the AIR IN (REGULAR) mode

- Aim the water flow at the space between the teeth and gum.
- Slowly shift the water flow along the gum.
- Clean the front and back of the teeth and behind the molars as well.
- Clean the areas of concern in the gum in the AIR IN (SOFT) mode.
- * Do not clean the periodontal pocket in JET mode, since the strength of water flow is too strong.

► Gum care

In the AIR IN (REGULAR) mode

- Stimulate the gum.
- Gently stimulate the areas of concern in the gum in the AIR IN (SOFT) mode.

Orthodontic nozzle

When cleaning crooked teeth and orthodontic bracket area with the brush and water flow

Lightly place the brush tip on the teeth and wiggle the brush.

► For crooked teeth

- Place the brush tip on the area that is difficult for a toothbrush to reach and tends to remain unpolished.
- Wiggle the brush mainly in interdental spaces.



► For orthodontic bracket area

- Insert the brush to the gap between the teeth and the orthodontic bracket.
- Wiggle the brush tip.

After use

1. Open the water tank cap and discard the water.
 - Be sure to discard the water in the water tank every time.
2. Turn on the power switch and drain out the water.
 - Operate until no more water appears and then turn off the power switch.

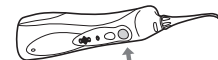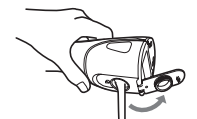
3. Stand the main unit in the charger.
4. While pressing the nozzle release button, remove the nozzle from the main unit and place it on the nozzle stand.

Cleaning

1. Pull out the water tank downwards.

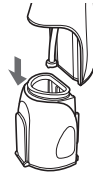
2. Clean with water.

- Remove the identification ring and clean with water.
- Do not wash with hot water exceeding 50 °C (122 °F).

- If using a detergent, use one that is suitable for dishwashing. Rinse well with water to ensure no detergent remains after cleaning.
- Wipe away moisture with a towel and leave to dry in a well-ventilated area.

Scrub the filter at the tip of the suction hose lightly with a soft brush.

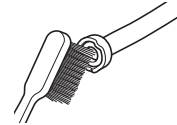
3. Wipe away the water with a towel etc. and dry.
4. Align and slide the water tank along the groove and install it securely.

Use a cloth to wipe away the stains on the charger.

- Never use thinners, benzine, or alcohol since it may cause malfunction or cracking/discoloration of the parts.
- Do not clean the charger with water.
- About once every six months wipe the prongs of the power plug with a dry cloth.

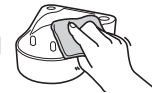Frequently asked questions

English

Question	Answer
Will the battery life be reduced if not used for a long period?	▶ When the appliance is not used for 6 months or more, the battery will weaken (leak battery fluid, etc.). Fully charge the battery once every 6 months.
Can the battery be charged before every use?	▶ Yes, but it is recommended that you charge the battery when it is empty. The life of the battery will vary depending upon usage and storage conditions.

Troubleshooting

Problem	Possible cause	Action
The appliance does not operate.	▶ You have just purchased the product or it has not been used for more than 6 months.	▶ Charge the unit for at least 23 hours. (See page 10.)
	▶ The charging time is short.	
The appliance can be used only for several minutes, even if it is charged.	▶ The main unit is not standing upright on the charger. (See page 9.)	▶ Stand the main unit upright so that its base contacts the charging section of the charger.
	▶ Battery life is over (approx. 3 years).	▶ The battery has reached the end of its life. (See page 15.)
Water pressure is low.	▶ The nozzle is clogged.	▶ Replace the nozzle.
	▶ Mouth of the nozzle is squashed.	
	▶ Filter is clogged.	▶ Clean filter. (See page 13.)
Water does not come out.	▶ The main unit is tilted too much while in use.	▶ Use with the main unit in a vertical state. (See page 11.)

If the problems still cannot be solved, contact the store where you purchased the unit or a service center authorized by Panasonic for repair.

Battery life

The battery life is about 3 years when charged about once a week (when used by one person).
If the operating time is significantly shorter even after a full charge, the battery has reached the end of its life. (The operating life of the battery may differ depending on usage or storage conditions)
Have the battery replaced by an authorized service center.

Removing the built-in rechargeable battery

ATTENTION:

A nickel metal hydride battery that is recyclable powers the product you have purchased. Please call 1-800-8-BATTERY for information on how to recycle this battery.


Remove the built-in rechargeable battery before disposing the appliance.

The battery is to be disposed of safely.
Please make sure that the battery is disposed of at an officially designated location if there is one.

This figure must only be used when disposing the appliance, and must not be used to repair it. If you dismantle the appliance yourself, it will no longer be waterproof, which may cause it to malfunction.

- Remove the appliance from the charger when removing the battery.
- Press the power switch to turn on the power and then keep the power on until the battery is completely discharged.
- Perform steps ① to ⑥ and lift the battery, and then remove it.
- Take care not to short circuit the positive and negative terminals of the removed battery, and insulate the terminals by applying tape to them.

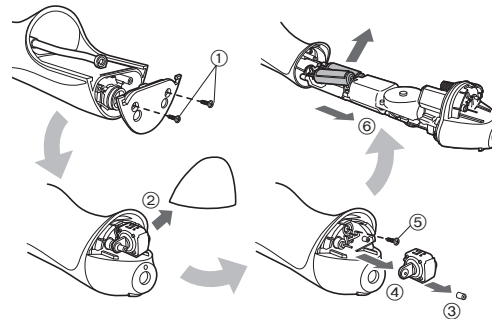Specifications

Power source	See the name plate on the product (Automatic voltage conversion)
Charging time	Approx. 15 hours
Power consumption	Approx. 1 W
Operating time	Approx. 15 minutes* (at 20 °C (68 °F) when fully charged)
Tank volume	Approx. 130 mL (Nozzle (for general use): approx. 35 seconds*) (Orthodontic nozzle: approx. 30 seconds*)
Dimensions	Main unit: 197 (H) x 59 (W) x 75 (D) mm (7.75 (H) x 2.32 (W) x 2.95 (D) in) (Not including nozzle) Charger: 40 (H) x 93 (W) x 98 (D) mm (1.57 (H) x 3.66 (W) x 3.85 (D) in)
Mass	Main unit: Approx. 305 g (Including nozzle) Charger: Approx. 180 g

* When using AIR IN (REGULAR) mode
This product is intended for household use only.

Changing the nozzle (Sold separately)


- Optional nozzles with identification ring to allow several persons to share the main unit.
- **Nozzle (for general use)**
- If the tip of nozzle (for general use) gets deformed, please replace with a new nozzle.
- The nozzle is a consumable. The nozzle (for general use) should be replaced with a new one every 6 months or so due to hygiene reasons, even if its shape is unchanged.

New nozzle (tip) Deformed (worn) nozzle (tip)

► Orthodontic nozzle

- Replace the orthodontic nozzle in approx. 3 months. If the brush bristles are bent out, replace it even within 3 months.

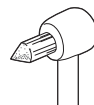New nozzle

Nozzle with bent out bristles

Part number	Replacement Nozzle for EW1213	
EW0950	Nozzle (for general use)	
WEW0984	Orthodontic nozzle	

Federal Communication Commission Interference Statement

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules and the limits for ISM equipment, pursuant to Part 18 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

FCC Caution: To assure continued compliance, install and use in accordance with provided instructions. Any changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate this equipment.

This device complies with Part 15 and Part 18 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Personal Care & Beauty Care Products

FOR USA ONLY

30-Day Quality Satisfaction Guarantee

If you are dissatisfied with any Panasonic Personal Care & Beauty Care Products for any reason, simply return it to the place of purchase with a dated proof of purchase, in the original packaging, with all accessories, parts and instructions, within 30 days of the date of purchase, for a full refund, or call Panasonic at 1-800-338-0552. Abuse or misapplication of this product voids the guarantee.

PRODUCT REGISTRATION

(U.S. customers only)

Go to www.panasonic.com/register and register your product now to receive these benefits:

Safety Notification

We contact our customers in case of recall

Owner Verification

Proof of purchase in case of product theft or loss

Efficient Warranty Service

In case there is a problem with your product

Special Offers and Promotions from Panasonic

**FOR YOUR CONVENIENCE, WE HAVE
ESTABLISHED A TOLL-FREE (IN USA ONLY)
NATIONAL ACCESSORY PURCHASE HOT LINE AT
1-800-332-5368.**

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© Panasonic Corporation 2020
Panasonic Corporation of North America
Panasonic Manufacturing (Thailand) Co.,Ltd.